Safeguarding hygine where it matters most





kikkoman<sup>®</sup>

#0001

2260







Lumitester Smart

Lumitester™

Smart



Hygiene monitoring in three easy steps



Visualize the data with App



Access from anywhere with cloud storage



With A3 Detection

ATP + ADP + AMP

More Reliable and Sensitive A3 Finds What Others Miss



www.hyserve.com



info@hyserve.com



+49 (0) 8846 - 13 44





With the App, anyone can check cleanliness easily



Easy-to-use, results within 10 seconds



## Swab

Swab the sample and insert the swab into main body



## Measure

Shake well and insert into the Lumitester





Analyze Measurements (RLU) displayed within 10 seconds



App to continuously monitor multiple data

Recorded data can be displayed in time-series by registering with our specialized App. Inspection pass rates are automatically graphed and improvements can be visualized. Employees will be more conscious of hygiene, and will maintain a high standard of environmental sanitation, which will help build greater trust in your business



Display time-series data for each inspection




Visualize overall inspection scores through graphs



Centralized multi-site data with cloud data management

Data will be stored in the cloud and can be accessed from anywhere in the world. This will allow fast detection and response to any problems that may arise.

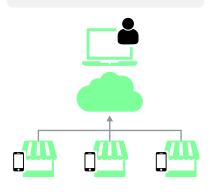

**Lumitester** ™ Smart Product code: 1902653

Measurement time 10 seconds

Data output RLU (Relative Light Unit)

Power 2 alkaline or nickel hydride rechargeable batteries (AA)

Included Accessories 2 alkaline batteries (AA), Cleaning brush, USB cable, Strap, Quick Manual

- Use for hygiene monitoring and other applications where detection of organic residues needs to be performed.
- Lucipac A3 and Lumitester Smart should not be used for counting general living bacteria or detecting specific pathogens.

LuciPac™ A3 Sui		Product code: 1702671 Product code: 1702672	
Storage Condit	ion	2-8°C (Do not freeze)	
25°C		14 days (Unopened)	
30°C		5 days (Unopened)	
Expiry		15 months after manufacturing date	



HyServe GmbH & Co KG. Hechenrainerstrasse 12 Tel. 0049 – 8846-1344

info@hyserve.com

D- 82449 Uffing am Staffelsee Fax 0049 – 8846-1342 www.hyserve.com